


Retrieval Bag

For easy collection and flexible extraction of tissue specimens



PROVIDING SPECIALISTS WITH
INNOVATIVE
HEALTHCARE SOLUTIONS

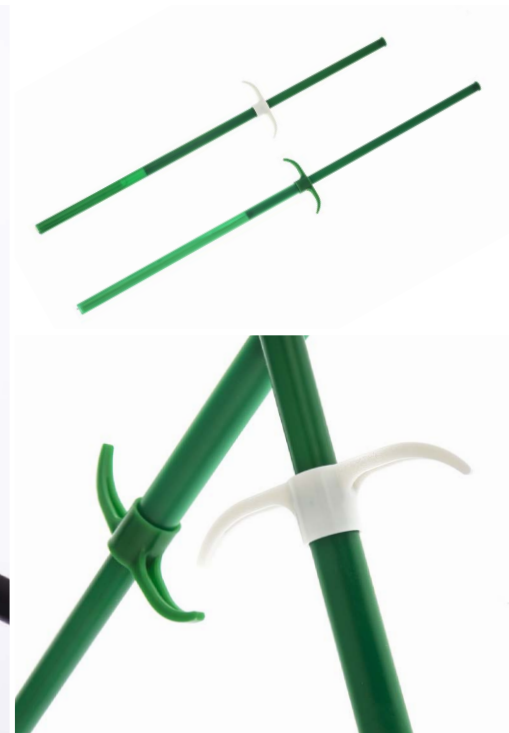




WITH THE CENTRAL FOCUS ON
PATIENTS' SAFETY AND USERS' COMFORT



Retrieval Bag

The Retrieval Bag is a disposable device used as a receptacle for the collection and extraction of tissue specimens, such as the appendix, gallbladder, ovaries, fibroid tumors, other tissues and calculi during laparoscopic surgical procedures. The Retrieval Bag is comprised of a flexible plastic bag slipped on the shape memory wire, which facilitates opening after deployment in the body cavity and maintaining in a fully open position. The bag is located in a plastic deployment tube but not permanently attached to it that does not block the trocar and allows to use it after the bag deployment. The bag can stay inside the body cavity and be filled with tissue a few times during surgical procedure and removed just before surgery is completed.



REF	COLOUR	VOLUME	TROCAR DIAMETER	PCS PER BOX
0208-RBM200		200 ml	10 mm	30
0208-RBM400		400 ml	10 mm	30
0208-RBM800		800 ml	10 mm	30
0208-RBM1200		1200 ml	10 mm	30
0208-RBM1500		1500 ml	10 mm	30

 **Grena Ltd.**
1000 Great West Road,
Brentford, London TW8 9HH,
United Kingdom

LONDON CHAMBER
COMMERCE *of* AND INDUSTRY

CUSTOMER SERVICE CENTER
NOTTINGHAM OFFICE, UK
☎ + 44 115 9704 800
☎ + 44 115 9704 800
✉ export@grena.co.uk
www.grena.co.uk

