

# Blue Endo-Hook Electrodes


Monopolar, Single Use, for Laparoscopic Procedures

*“For safe and fast performance during laparoscopic procedures”*

- Available with various tip versions and lengths
- Fits all diathermy handles made for 2.35mmØ electrodes
- Insulated to tip for better safety
- Hexagon connector for firm fitting
- Cost effective




**Blue Surgical**  
your OR partner

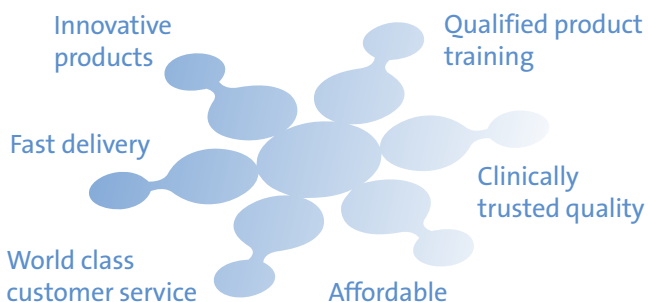

## Blue Endo-Hook Electrodes

Blue Endo-Hook Electrodes are designed for comfort of use in Laparoscopic procedures. Attach to any diathermy handle made for 2.35mmØ electrodes. Simple and convenient. The lightweight 5mm electrode shaft is fully insulated, all the way to the tip, ensuring optimum safety for patient and surgeon.

Available with standard L-tip, which is used in most Laparoscopic procedures, or a larger L-tip with flat top specially designed for Thorax procedures. L-tip electrode is also available attached firmly on a diathermy handle from the start, in either 33cm or 45cm length. Ready to use.

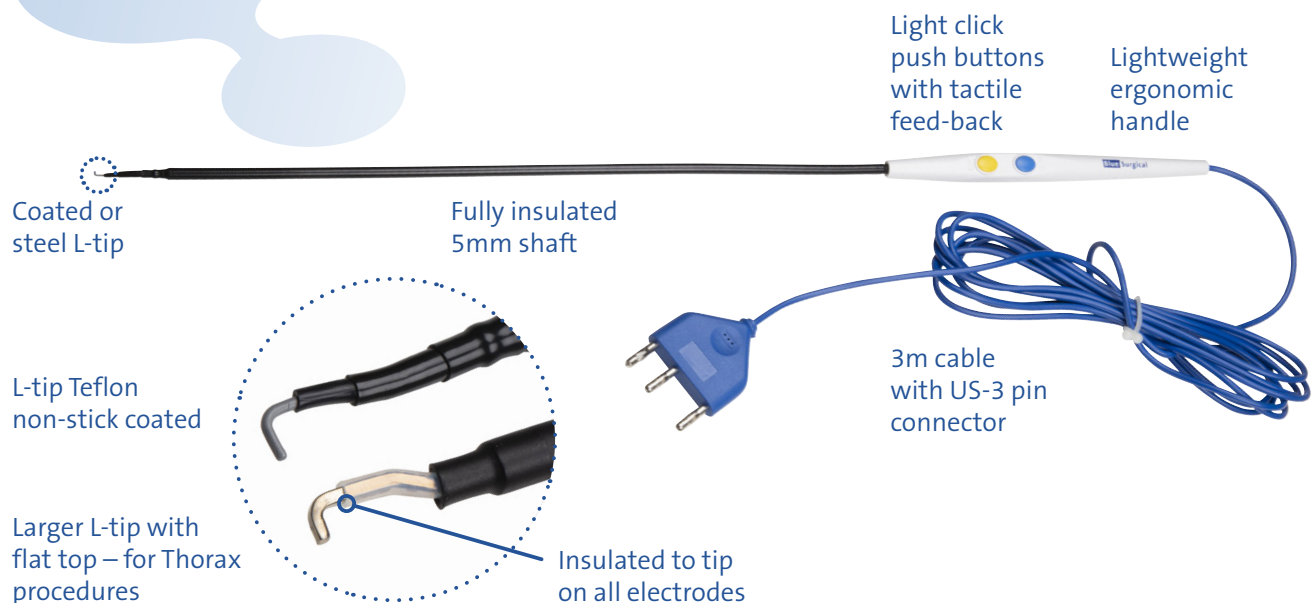

REF	Description	Overall length	Pcs / box
<b>Electrodes alone (without diathermy handle):</b>			
B-46500	Blue Endo-Hook Electrode, L-tip Teflon coated, single use	5mm x 33cm	10
B-46600	Blue Endo-Hook Electrode, L-tip Teflon coated, single use	5mm x 45cm	10
SW133407214S-BM	Blue Endo-Hook Electrode, large L-tip, with flat top, steel, single use for Thorax procedures	5mm x 30cm	12
<b>Diathermy handle alone:</b>			
B-40400	Blue Diathermy Handle, finger-controlled with 70mm steel electrode, 3m cable, US-3 pin connector, single use		25
<b>Complete ready to use set: Diathermy handle + Endo-Hook Electrode</b>			
B-46800	Blue Endo-Hook, finger-controlled L-tip Teflon coated, with 3m cable and US 3-pin connector, single use	5mm x 33cm	10
B-46810	Blue Endo-Hook, finger-controlled L-tip Teflon coated, with 3m cable and US 3-pin connector, single use	5mm x 45cm	10

**Blue Surgical**  
your OR partner

Exclusive International Distribution by:

Blue Surgical ApS • International Sales Division • Virumvej 64 • DK-2830 Virum • Denmark  
Phone: +45 21 19 22 30 • E-mail: info@bluesurgical.com • www.bluesurgical.com