


# ChorionSet

Transabdominal Chorionic Villus sampling set  
for fetopathies determination

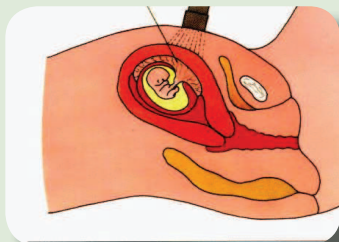

**ChorionSet (PBUT)** is a set containing two devices:  
**A- coaxial needle**  
**B- needle for villus sampling**

The coaxial needle allows the sampling needle to reach the sampling site through the abdominal wall.  
The two needles are perfectly coupled and allow an optimal positioning.

The sampling needle has a cannula with special sharpening called **“UltraPoint”**, specifically designed to enable a maximum penetration capacity, keeping safe the fetus.

The **“echo” surface treatment** allows a correct placement of the needle and its constant visualization under ultrasound guidance.

The needle is coloured according to the **interantional colour codes** for easy gauge size identification.

The introducer needle is provided with **wings** to allow a better grip.

Code	Gauge	Outer diameter	Inner diameter	Length
PBUT2015	18G (introducer)	1.27 mm	1.05 mm	90 mm
	20G (needle)	0.90 mm	0.70 mm	150 mm
PBUT2020	18G (introducer)	1.27 mm	1.05 mm	140 mm
	20G (needle)	0.90 mm	0.70 mm	200 mm
PBUT2022	18G (introducer)	1.27 mm	1.05 mm	160 mm
	20G (needle)	0.90 mm	0.70 mm	220 mm
PBUT2222	19G (introducer)	0.90 mm	0.70 mm	160 mm
	22G (needle)	0.70 mm	0.50 mm	220 mm



Optionally NeedleVise® safety device can be added **(PBUTR)**

CHORIONIC VILLUS  
SAMPLING



M.D.L. Srl  
Via Tavani 1/a - 23014 - Delebio (So) Italy  
Tel +39 -0342 682130 Fax +39 0342 691316  
email : uc@mdlsrl.com  
web : www.mdlsrl.com

CE  
0426