


*The solution  
disposing fluids* + DSXr



**Drain Star  
Revolution**

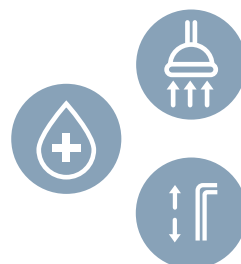

[www.drainstar.it](http://www.drainstar.it)

*Safe, innovative  
and eco sustainable  
unit for hazardous  
waste management.*

*Drain Star Revolution, the next generation  
in surgical fluid collection and disposal systems,  
is the result of over 15 years experience.  
The high volume capacity allows uninterrupted  
performance for OT staff, managing fluids  
by a portable unit totally closed.*



## + BENEFITS

Eliminate canister replacement  
& disposal

Reduce solid waste disposal costs

Operate with safe, continuous suction  
for any surgeries

Environmentally friendly solution  
to bioburden removal during surgery

Battery powered for use in ward as I.C.U.  
(CVH), dialysis or catheterized patients

# + FUNCTIONS


# DSXr



**Suction**  
Two independent suction lines low vacuum




**I.V. Pole**  
Lifter with electrical height adjustment

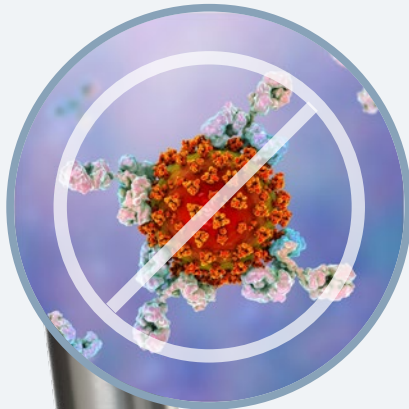

**Bacterial neutralization**  
of the collected fluids

Automatic disinfection cycle

High accuracy volume measurement

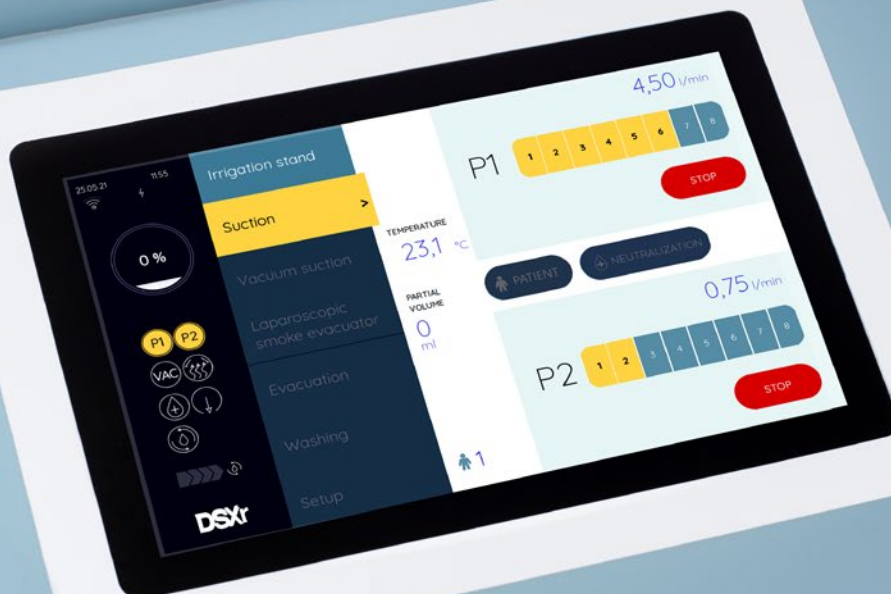

# DSXr



## TECHNICAL SHEET

Dimensions	cm 55 x 55 x 160 h
Weight	72 kg
Structure	Stainless steel
Wheels	125 mm diameter, two with brake
Supply	Voltage 110V – 230VAC 50/60 Hz Cyclical battery 12V 14Ah AGM
Capacity	62 litres
I.V. Pole	Maximum load capacity 30 Kg, regulating range 1.500 – 2.650 mm (Optional)
Low vacuum suction	Two independent suction line
Evacuation	20 lt/min
User interface	10.1" Display Touchscreen
Mode of operation	Continuous/Intermittent by mains supply or battery
Data memory	Up to 42 patients
Temperature	By display with configurable alarm
High level disinfectant reservoir	1,5 lt internal container
Decontamination cycle	Totally automatic



# Drain Star Revolution+

Drain Star by M.A.S. MEDICAL S.R.L.

**Registered Office**  
Corso Siracusa 57  
10136 - Torino - Italy

**Operating Office**  
Via Vittorio Alfieri 16  
12035 - Racconigi (CN) - Italy

**MAS medical**

Ph./Fax. +39 0172 84981  
info@drainstar.it  
www.drainstar.it